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
Method of Thawing	Time Required to Thaw	Thawed Cartridge Storage Life (2 – 8°C or on ice)	A cartridge should not be re-frozen after thawing
Room temperature water bath (opened)	~ 1 hour	24 hours ¹	
2 – 8°C (unopened)	8-12 hours	7 days	

Table 1. Guidance for reagent kit thawing and storage after thawing.

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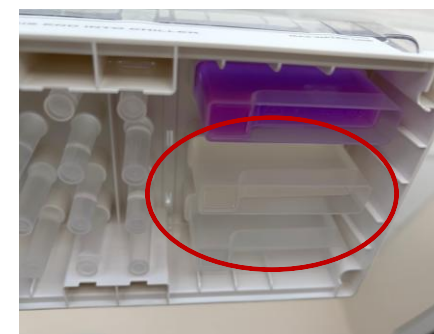Dilución y desnaturalización

Library Prep Kit	Initial Library Pool Concentration	Denatured Library Concentration After Addition of 990 µl of HT1
Nextera XT, NEB Ultra II FS, QIAseq FX, & KAPA HyperPlus	2 nM	10 pM
Illumina DNA Prep, PulseNet Rapid Prep, Nextera XT, NEB Ultra II FS, QIAseq FX, & KAPA HyperPlus	4 nM	20 pM
Illumina DNA PCR-Free	2 nM	20 pM

Table 2. Concentrations for the initial library pool and the denatured pool for PulseNet validated library preps

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Table 3. Recommended loading concentrations for different library prep kits. ¹Not determined.

Final Loading Concentration	Denatured Pool Concentration = 10 pM		Denatured Pool Concentration = 20 pM	
	Required Volume of HT1	Required Volume of Denatured Pool	Required Volume of HT1	Required Volume of Denatured Pool
8 pM	200 µl	800 µl	600 µl	400 µl
10 pM	NA	NA	500 µl	500 µl
12 pM	NA	NA	400 µl	600 µl
15 pM	NA	NA	250 µl	750 µl

Table 4. Required pool and HT1 volumes for dilution to desired final concentration.

Dilución y desnaturalización

PNL35.W1 Illumina DNA Prep Workbook

Calculate molarity for all instrument types:	
Pool concentration (Qubit value, ng/μl):	15
Molarity (nM) = (conc/(660g/mol x 1000bp))x 10 ⁶ :	22.73
Volume of pool needed (μl) = (nM)(x)=(4nM)(50μl):	8.80
Volume of RSB needed (μl) = 50 - volume of pool:	41.20
Result is 50 μl of 4 nM pooled libraries	

Calculate molarity for all instrument types:	
Pool concentration (Qubit value, ng/μl):	15
Molarity (nM) = (conc/(660g/mol=(PoolConcentration/(1000*660))*10 ⁶	
Volume of pool needed (μl) = (nM)(x)=(4nM)(50μl):	8.80
Volume of RSB needed (μl) = 50 - volume of pool:	41.20
Result is 50 μl of 4 nM pooled libraries	

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MiSeq run, check if yes	
5 μl of 4 nM pool + 5 μl of 0.2 N NaOH + 990 μl of HT1 = 20 pM of denaturated pool	
Desired Final Loading concentration (12 pM to 20 pM):	12
Volume of 20 pM library pool:	600
Volume of HT1 buffer:	400
Result is 1000 μl of library pool at the desired loading concentration	
If spiking PhiX: Dilute 10 nM PhiX to 2 nM: 1 μl of 10 nM + 4 μl of RSB	
Denature: 5 μl of 2 nM PhiX + 5 μl of 0.2N NaOH (now at 1 nM)	
Dilute to 20 pM: 490 μl of HT1 + 10 μl of 1 nM PhiX	
For v2 chemistry, dilute PhiX to 12.5 pM: 62.5 μl of 20 pM PhiX + 37.5 μl HT1	
Result is 100 μl of PhiX at desired loading concentration	
Percentage of PhiX spiked (if applicable):	

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